



ClorDiSys

"The Chlorine Dioxide People"

Providing you with gaseous chlorine dioxide solutions for your decontamination and sterilization needs

CEILING PLATE

Description:



The Ceiling Plate works in conjunction with the Distribution Plate (PLAT-DISTXX-001) to provide a method of decontaminating a single room from a remote location. The ceiling plate is mounted in the room and connected to the Distribution Plate via tubing and signal cables to communicate with the CD Gas Generator located next to the Distribution Plate at a central location.

The ceiling plate is to be securely screwed and sealed to the ceiling to prevent gas from leaking during the decontamination process.

The Ceiling Plate can be configured with a Pressure Relief valve for rooms where pressure build up may be an issue, such as BSL3-Ag facilities and other extremely tight rooms.



Ceiling Plate
(Part #: PLAT-CEIL/PR/RH-001)
Unit Dimensions: 13.0" x 13.0"
Unit Weight: 13 lbs

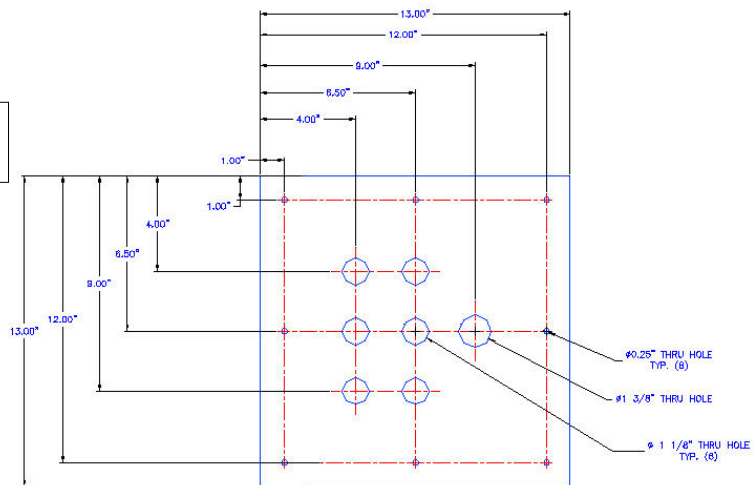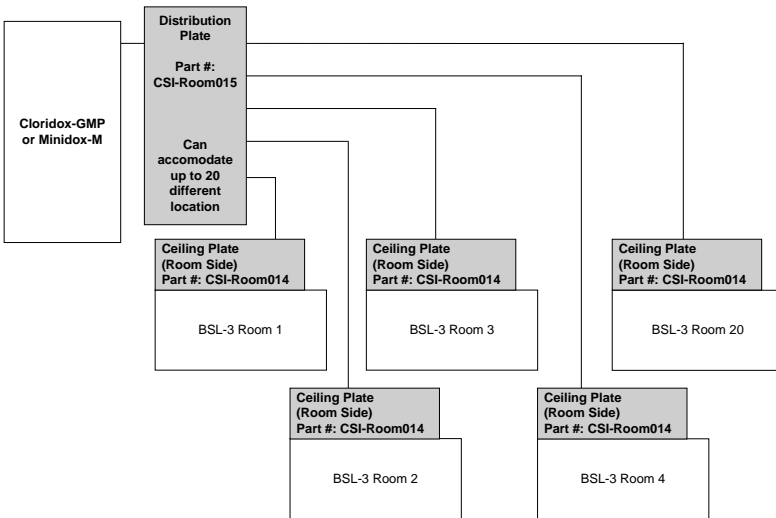

Product Number	Description
PLAT-CEIL/PR/RH-001	Door Plate, Pressure Relief and RH Probe included.
PLAT-CEIL/RH-001	Door Plate, RH Probe included (No Pressure Relief)
PLAT-CEIL-001	Door Plate, RH Probe not included and no Pressure Relief



P.O. Box 549, Lebanon, NJ 08833-0549 Tel: (908) 236-4100 Fax: (908) 236-2222

www.clordisys.com